


Q
U
A
L
I
T
Y



Also ask to see the
Entrée and Recumbent

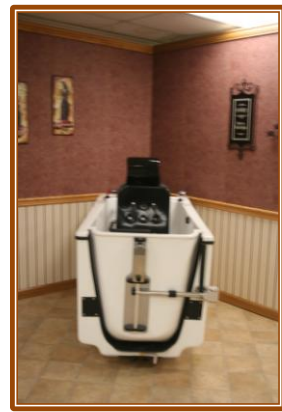

W
A
R
R
A
N
T
Y

3 Year
Standard

10 Year
Extended
Available

Standards & Certifications:

Our Mission Statement is: "Through Good Manufacturing Practices, **MasterCare** is committed to design, manufacture and distribute the highest quality bathing systems for physically impaired individuals throughout North America. Our systems are backed by professional customer and technical service.

MasterCare's bathing systems are tested and certified to meet or exceed the following standards:

- Electrical Standards:** **CAN/CSA C22.2 No. 218.2-93 (R2004)** Hydromassage Bathtub Appliances. **UL 1795** Hydromassage Bathtubs, Third Edition, Including revisions through to Sept. 22, 2006.
- Plumbing Standards:** **Warnock Hersey, CSA B45.5** (2008); **ASME A112.19.15** (2005); **CSA B45S1** (2004); **ANSI Z124.1.2** (2005).
- State Standards:** Please call 800-798-5867 for information on your state standards.
- FDA Establishment Registration No.** 3003333691



Complimentary tub kit, install inspection & in-service training provided by your local area certified representative—no hidden costs for accessories.



MasterCare

2071 14th Avenue
Columbus, NE 68601

800-798-5867 ph
402-563-9102 fax
www.mastercarebath.net



MasterCare



MB-80-R Front
Entry Bathing



Made in the USA



I
N
T
E
G
R
I
T
Y

FEATURES

MB-80-R Integrity Bath

Standard Features:



Pre-fill reservoir fills the bath quickly

Analog temperature readouts

Deck mounted anti-scald high volume Thermoscopic mixing valve

Entire plumbing assembly panel removes for easy servicing

Concealed shower hose storage when not in use

Two inch overflow

Front-entry system with optional multi-position reclinable chair.

Ergonomic design of the chair provides the deepest immersion, for this type of system, on the market today.

Headrest eliminates resident neck strain and assists bath aide in shampooing.



Flow through holes for increased water circulation.

Chair locks to the base for transport and to rail system for bathing. An individual only needs to be transferred once.

The adjustable leg rest lifts an individual's legs for transfer and provides easy access to lower extremities for the caregiver.



Optional Scale; weighs up to 550 lbs

Locking casters for safety and security

Easy access service doors also function as lockable liquids storage with shelves

Disinfecting system; does not require proprietary shaped bottles

Easy open drain handle

On-board liquids dispensing

Two inch pop-up drain



Create the perfect seal each and every time.

Simple as 1 - 2 - 3

- ❖ Rotate the door inside the bath.
- ❖ Push down to lock.
- ❖ Rotate cam arm up.



1



2



3



Casual water drain system to divert any minor leaks and residual water into drain.

Please refer to **MasterCare's** heated variable speed **BathAire^(SM)** brochure for more information on the **only closed-loop** air system on the market today.



Aromatherapy sachets are inserted into an aromatherapy canister that carries the scent into the BathAire^(SM) system.

GFCI Reset located behind service lockable storage door.

Installs



Westview Care Center
Britt, IA



Laurels of Bon Air
Richmond, VA



Cross Keys Village
Oxford, PA